


STAALAFNAME CEREBROSPINAAL VOCHT

voorbereiding	Deze punctie moet in steriele omstandigheden uitgevoerd worden. De voorbereiding gebeurt conform de verpleegkundige procedures ter zake.
staalafname	Afname van cerebrospinaal vocht dient door een ervaren arts uitgevoerd te worden! Neem steeds meerdere steriele, gegradueerde polypropyleen tubes af, en nummer deze in volgorde van afname.
aandachtspunten	Dit is gewoonlijk een kritisch en DRINGEND staal: Breng het staal ONMIDDELIJK naar het labo, en geef het af in handen van een labo medewerker. Transport via buizenpost wordt afgeraden.

NEEM HET JUISTE RECIPIËNT VOLGENS DE AANGEVRAAGDE ANALYSES VOOR STAALTRANSPORT.

- AANVRAAG ALGEMENE CHEMIE - CYTOLOGIE - BACTERIOLOGIE

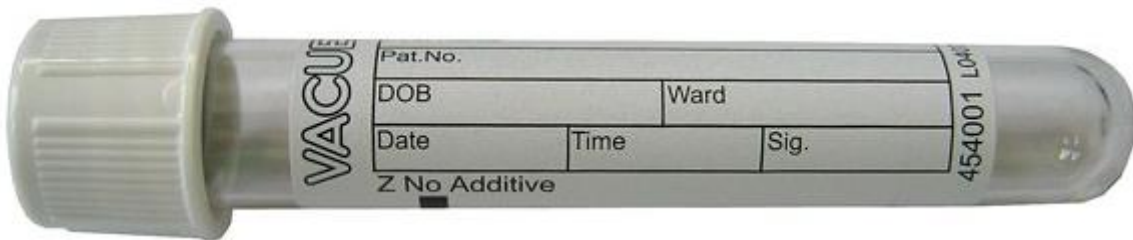

- Afname dient te gebeuren in onderstaande microtubes bij specifieke aanvraag CSV biomerkers (Ziekte van Alzheimer, Creutzfeldt-Jakob)

Lumbar puncture



- Lumbar puncture (LP) at level L3/L4 or L4/L5
- Discard the first 2ml of CSF


Collection



- Collect at least 2x 1ml-1,5ml CSF directly into 1,5ml polypropylene (PP) Low Bind (LoB) tubes*

*Order nr 72.703.600 - 1,5ml micro tube SARSTEDT AG & Co. KG

Storage / Transport

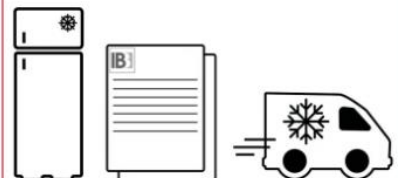

- Ship fresh CSF within 7 days (max 4°C) or 48h (max 22°C)
- Frozen CSF must be kept in frozen condition while shipping
- Ship plasma always frozen, preferably at -80°C
- Include 2 forms: our Biomarker request sheet and your lab subcontracting form.